

STERISAMPLE®



...pure sampling

Sample port

The hypodermic sample port makes it possible to take samples of the product from a vessel in an aseptic manner without disrupting production. It is suitable for a wide range of processes/production lines within the pharmaceutical sector. Prior to sampling, the port is sterilized with alcohol and then sampling takes place by simply penetrating the special silicone membrane with a needle.

Features

Easy to use. Easy to clean/sterilize. Flanged on to a standard TC-flange or the Steriflange. Adaptation to a syringe or to a vessel by flexible hose depending on the sample volume. Cost effective sampling system by using standard syringes. Autoclavable.

Technical data

Product contact surfaces are finished to $Ra < 0.5 \mu m$ or better.

Silicone diaphragm complies with FDA regulations.


Traceability-certificates according to EN 10204 (3.1B) are always supplied with Sterisample.

Manufactured of AISI 316L stainless steel.


Applications

General Sampling, Cell Culture Sampling and High Purity sampling.

Ways of taking samples from vessel wall

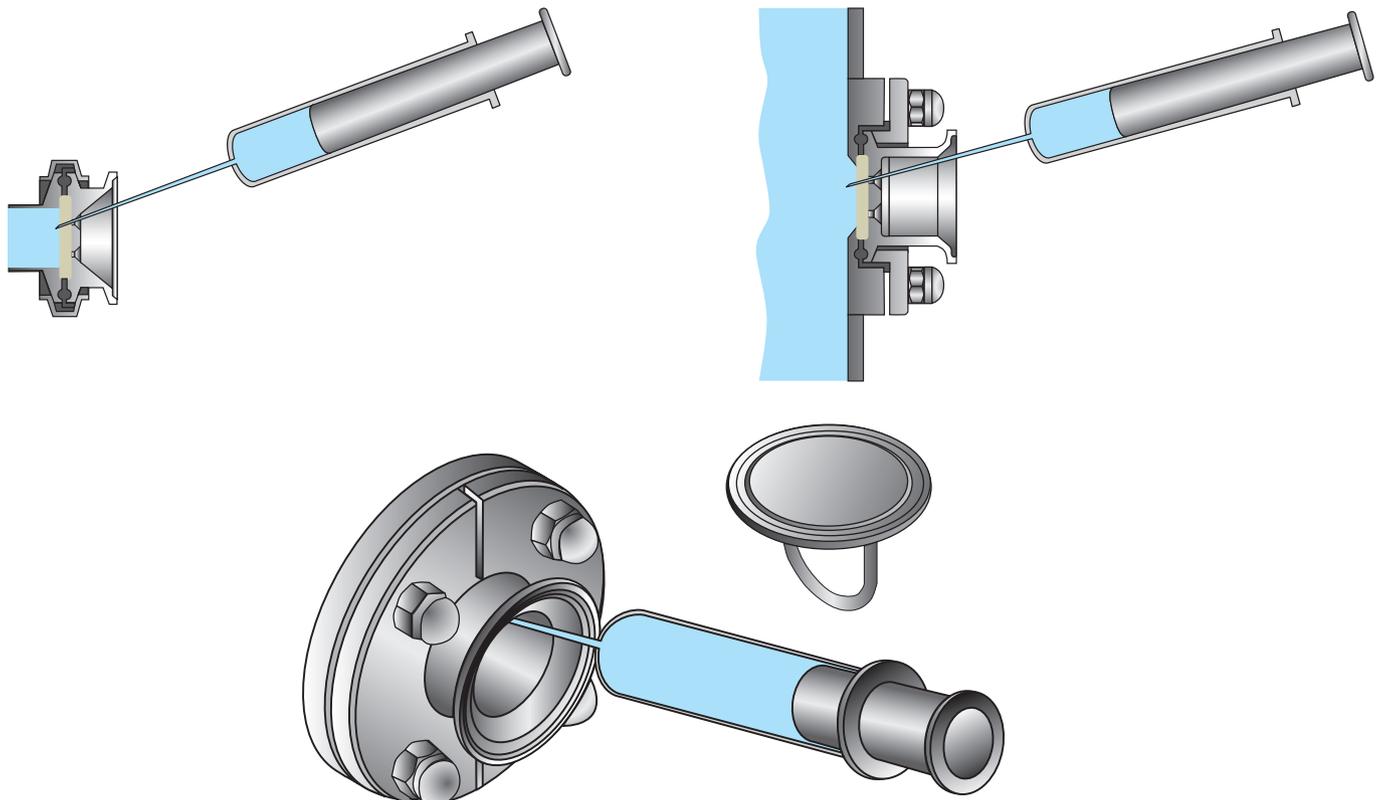


Sterisample short, for standard TC-flange Connection

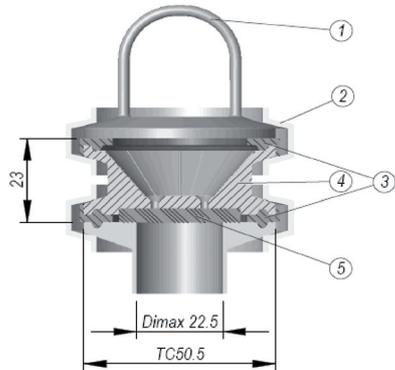

110336

6	1	Nut for TC clamp	110339
5	1	Membrane	110034
4	1	Sample port 50.5	110332
3	2	TC-gasket 50.5/38	110334
2	2	TC-clamp 50.5/38	110014
1	1	Blind endcap 50.5/38	110335
Item	Qty	Denomination	Part No

Sterisample long, for Steriflange Connection

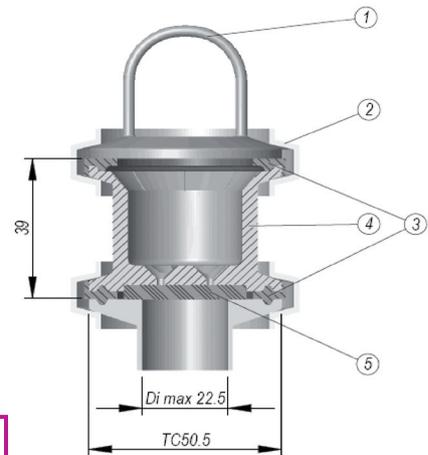


110337


6	1	Nut for TC clamp	110339
5	1	Membrane	110034
4	1	Sample port 50.5	110333
3	2	TC-gasket 50.5/38	110334
2	1	TC-clamp 50.5/38	110014
1	1	Blind endcap 50.5/38	110335
Item	Qty	Denomination	Part No

Recommended Steriflange for use with Sterisample


SF ASTM M-100/12

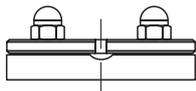


SF ASTM M-100/25

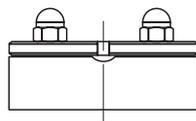

SF ASTM M-100/38

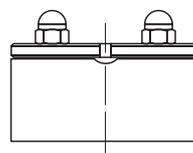

SF ASTM M-100/51

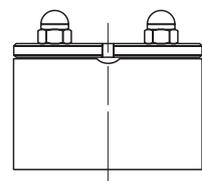

Sterisample fits in the following clamp connection Steriflange ASTM M-100 (metric), and ASTM I-100 (inch), and can be supplied in the following four sizes, 12, 25, 38, and 51 mm.



STERIDOSE SALES

www.steridose.com

